

ANKALAPS
SAFE BAG
Specimen Retrieval Bag



H Series




HW Series



W Series



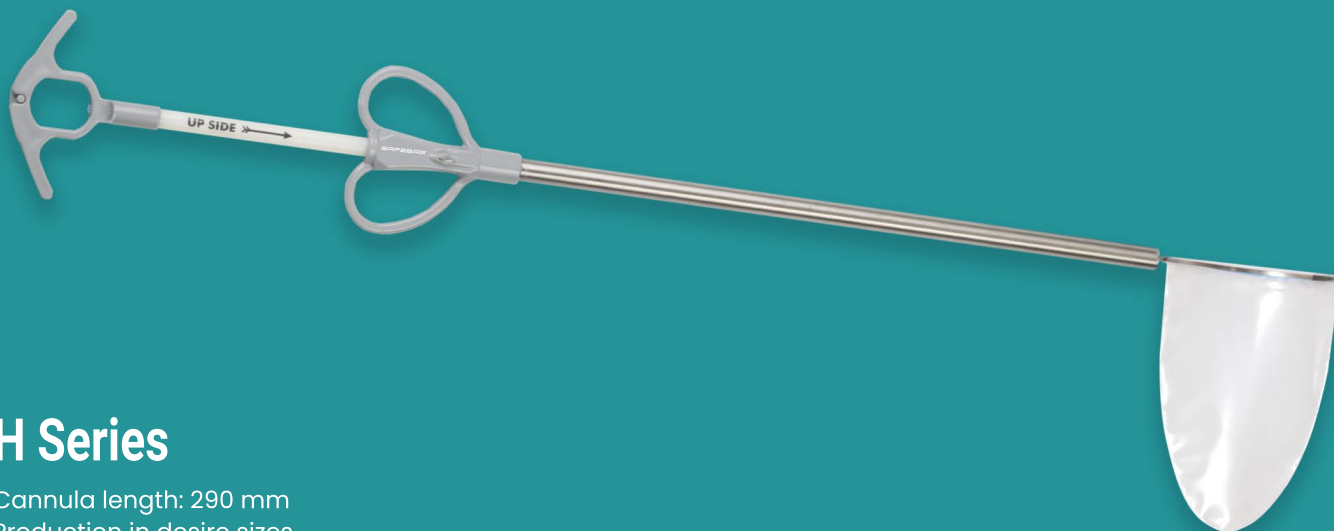

H Series

Cannula length: 290 mm
Production in desire sizes

Features

Endobags of H Series can be closed off more than one in abdomen. They are produced two different design as structure of asymmetric mouth and circular mouth. Endobags are manufactured from polyuretan that stretches 80 and 100 micron.

Extra wide and asymmetric mouth structure facilitates insertion of specimen

Intuitional design facilitates the extraction of specimen

Extraordinary durability and closure

Maximize the specimen's safety
Minimize the risk of contamination
Provides complete isolation for cancerous tissues

Antibacterial, Non-toxic and Non-allergic

Provides safe us efor all patient groups

Flexible and Transparent

Visible specimen

H Series Dimension Table

Model No	Size	Dimension	Pieces/Box	Volume
S10H	Small	65 mm x 150 mm	20 Pcs	200 ml
M10H	Medium	75 mm x 200 mm	20 Pcs	400 ml
L10H	Large 10	100 mm x 200 mm	20 Pcs	800 ml
L12H	Large 12	125 mm x 230 mm	20 Pcs	1500 ml
XL15H	X-Large	130 mm x 250 mm	20 Pcs	2000 ml



HW Series

Cannula length: 290 mm
Production in desire sizes

Features

Endobags of HW Series open automatically by aid of shape memory metal strip, located at the mouth of bag. They can be removed from metal hoop by pulling back the ring, connected to the inner shaft. Hereby, endobags do not engage with trocar ports during operations.

Extra wide and asymmetric

mouth structure facilitates insertion of specimen

Easy and Effective Usability

Bag can be easily separated from the metal hoop by pulling back the inner shaft.

The mouth of the bag can be shrink by pulling back the fused ring rope, located at the end of the inner shaft.

Does not leave particle during the separation of metal hoop.

The rope can be cut with hidden blade, located at the ring of inner shaft.

Allows the use of other tools without blocking the trocar.

Extraordinary Durability and Transparent Bag Structure

Protects both surgeon and patient from invasive reaction taken from specimens.

Prevents the risk of transmission of malignant cells.

Large specimens can be cut into small pieces.

HW Series Dimension Table

Model No	Size	Dimension	Pieces/Box	Volume
S10HW	Small	65 mm x 150 mm	20 Pcs	200 ml
M10HW	Medium	75 mm x 200 mm	20 Pcs	400 ml
L10HW	Large 10	100 mm x 200 mm	20 Pcs	800 ml
L12HW	Large 12	125 mm x 230 mm	20 Pcs	1500 ml
XL15HW	X-Large	130 mm x 250 mm	20 Pcs	2000 ml



W Series

Cannula length: 290 mm
Production in desire sizes

Features

Endobags of W Series are design for simple, memory and safe minimally surgical.

- ✔ Protects both surgeon and patient from invasive reaction taken from specimen.
- ✔ Through chrome-nikel multifilament, located at the mouth of bag, endobag is opened automatically and ready to use. After putting the specimen, the mouth of the bag closed easily with the help of a holder
- ✔ Can be opened and closed many times
- ✔ Allows the use of other tools without blocking trocar
- ✔ Large specimens can be cut into smaller pieces inside the endobag. Prevents the spreading risk of malignant cells. Produced from transparent, strong, non-allergic antibacterial polyuretan film. Available for use inside the protective cannula.

W Series Dimension Table

Model No	Size	Dimension	Pieces/Box	Volume
S10W	Small	65 mm x 150 mm	20 Pcs	200 ml
M10W	Medium	75 mm x 200 mm	20 Pcs	400 ml
L10W	Large 10	100 mm x 200 mm	20 Pcs	800 ml
L12W	Large 12	125 mm x 230 mm	20 Pcs	1500 ml
XL15W	X-Large	130 mm x 250 mm	20 Pcs	2000 ml